



## Next Generation Fluorescence Imaging

Smart Sensor Solutions

### General Information


This product is non-toxic, non-contagious and is not intended for human use. The plasmid DNA in the package is synthetic and only for research and development purposes. It does not present any danger for humans, animals or the environment. The product is only for in vitro studies and is not for commercial use.

### mtG-geNOP vector

This vector has not been completely sequenced. All provided information regarding the vector composition was compiled using the information from published literature, other sources together with partial sequences obtained by NGFI.

### Vector description

mtG-geNOP is a green fluorescent-based probe for imaging nitric oxide (NO•) within the mitochondrial matrix of intact mammalian cells. In order to express mtG-geNOP in cells of interest, 20 µg of purified endotoxin-free plasmid DNA coding for mtG-geNOP is provided. The plasmid coding for mtG-geNOP represents a mammalian expression vector with a strong viral promoter. For plasmid amplification in E.coli ampicillin should be used. 1 – 1.5 µg DNA is required for cell transfection in a single well of a standard 6-well dish following standard transfection procedures. Usually cells express high amounts of mtG-geNOP 24 – 48 hours after cell transfection. Standard optical filters for GFP imaging should be used. The vector can be also used as a source of mtG-geNOP coding sequence.

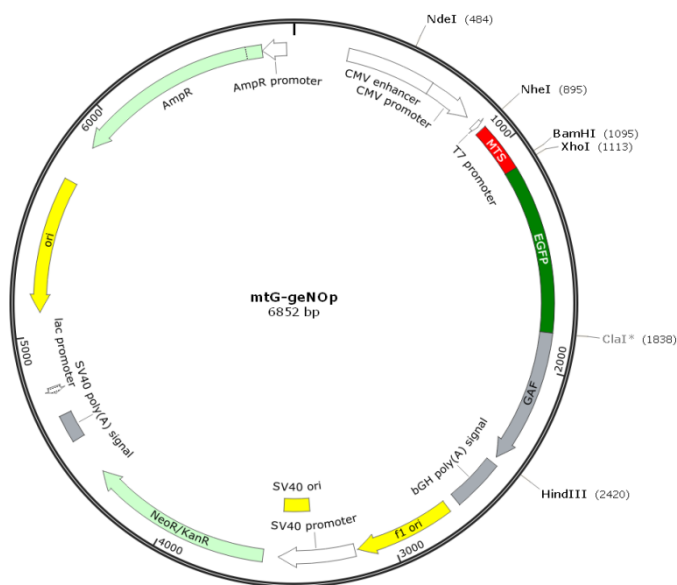

### Expression in mammalian cells

mtG-geNOP vector can be transfected into mammalian cells by any known transfection method. CMV promoter provides strong, constitutive mtC-geNOP expression in eukaryotic cells.

### Propagation in E. coli

Suitable host strains for propagation in E. coli include DH5alpha, HB101, XL1-Blue, and other general purpose strains. The vector confers resistance to ampicillin (100 µg/ml) to E. coli hosts.

NOTE: In order to supply the NO•-binding domain of mtC-geNOP with sufficient iron(II) cells need to be treated with the non-toxic iron(II) loading buffer (<http://www.ngfi.eu/product/ironii-booster-solution/>) for at least 10 minutes before imaging experiments.

### References

Eroglu E. (2016) "Development of novel FP-based probes for live-cell imaging of nitric oxide dynamics" (<http://www.nature.com/ncomms/2016/160204/ncomms10623/pdf/ncomms10623.pdf>)